




## GAIN INDEPENDENCE WITH THE MTG EZ-GRIPPER®

For decades, tetraplegics have been incapable of intermittent self-catheterization without help from a caregiver, so MTG set a goal of developing a closed system catheter for individuals with limited finger dexterity.

The MTG EZ-Gripper® is a sterile, single-use closed system catheter fitted with access finger holes for easier opening and handling and a patented gripping handle affixed to the outside of the collection bag.

Available in 12-16 FR sizes and 16" and 18" lengths.

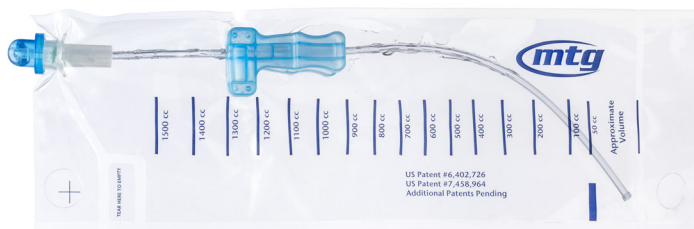

# MTG EZ-Gripper®

## No-Touch Closed System Catheter

### FEATURES & BENEFITS

#### Designed for Limited Dexterity

Our patented gripping handle attached to the collection bag allows for easy advancement, making it easier for individuals with limited hand dexterity to self-catheterize.

#### Features our EZ-Advancer® Valve

Our patented EZ-Advancer® valve secures the catheter tube after each forward advancement, preventing retraction into the bag during catheter insertion.

#### Designed to Reduce Risk of Recurrent UTIs

The no-touch closed system catheter and pre-lubricated introducer tip are designed to reduce the risk of recurrent urinary tract infections.

#### Collection Bag Protects Sterility

The self-contained 1500ml collection bag protects the sterility of the catheter tube and collects urine while measuring output. A sterile sample can be easily collected from the catheter port.

#### Made with Medical Grade Quality Materials

All MTG catheters use medical grade quality materials and are not made with natural rubber latex, BPA or DEHP.



MTG catheter kits include everything you need to self-cath. Full kits contain one sterile MTG EZ-Gripper® closed system catheter, pair of latex -free gloves, underpad, gauze pad, Benzalkonium chloride (BZK) clear, odor-free skin cleansing prep wipe and privacy bag.

## MTG EZ-Gripper® No-Touch Closed System Catheters

### Single Closed System Catheters

Item #	FR Size	Length	Description	Antiseptic Wipe
50112	12	16"	MTG EZ-Gripper® No-Touch Closed System Catheter, 12 Fr, 16"	NA
50114	14	16"	MTG EZ-Gripper® No-Touch Closed System Catheter, 14 Fr, 16"	NA

### Closed System Catheter Kits - Standard Size

Item #	FR Size	Length	Description	Antiseptic Wipe
52112	12	16"	MTG EZ-Gripper® No-Touch Closed System Catheter Kit, 12 Fr, 16"	BZK Wipe
52812	12	18"	MTG EZ-Gripper® No-Touch Closed System Catheter Kit, 12 Fr, 18"	BZK Wipe
52114	14	16"	MTG EZ-Gripper® No-Touch Closed System Catheter Kit, 14 Fr, 16"	BZK Wipe
52116	16	16"	MTG EZ-Gripper® No-Touch Closed System Catheter Kit, 16 Fr, 16"	BZK Wipe

HCPCS: A4353, Packaging: 50/Case

## PAIRS PERFECTLY WITH THE MTG EAGLE® BOARD

The MTG EZ-Gripper® closed system catheter can be paired with our MTG Eagle® Board, an adaptive reusable catheterization device, providing men with very limited dexterity even more independence to catheterize on their own.

- ✓ Makes One-Hand Catheterization Possible
- ✓ Penile Sizing Guide with Foam Pads for Precise Fit
- ✓ Thumb Holes for Easy Accessibility
- ✓ Washable and Reusable



**mtg Eagle Board**